

- ▶ Preamble

[1] 2 cells hidden

► **Data**


▶ 9 cells hidden

- Diffusion Coefficient

```

1 #data_list = list(zip())
2
3 fig, (ax1, ax2) = plt.subplots(2, figsize = (12, 8), sharex=False, sharey=True, layout='constrained')
4
5 ## List of tracks sorted by size. Length of list determines the readability of the plots
6 ##### Slice the list length by changing the integer i; here [0:i] and here [0:i]
7 zipped_list = zip(np.arange(tracks_by_size.shape[0], dtype='int')[0:5], tracks_by_size.astype('int')[0:5])
8
9
10 # plotting loop
11 for i,j in zipped_list: # np.arange(tracks_by_size.shape[0]) # np.arange(0, 5)
12     t = np.arange(1, data_links['DISPLACEMENT'].loc[data_links['TRACK_ID'] == tracks_by_size[i]].shape[0]+1)*0.05
13     d = np.array(data_links['DISPLACEMENT'].loc[data_links['TRACK_ID'] == tracks_by_size[i]]*6.25*1e-4
14
15     ax1.plot(t, d, label=j)
16     ax2.stackplot(t, d)
17
18 # combining all plots from the loop
19 ax1.legend(loc='upper center', bbox_to_anchor=(1.025, 1.1), ncol=1, fancybox=True, shadow=True, title='Tracks (ID)')
20 fig.supxlabel(r'Track length, $t$ [s]');
21 fig.supylabel(r'Displacement, $d$ [cm]');
22 fig.suptitle(r'Particle Displacement Versus Time')
23 plt.ticklabel_format(style='scientific', axis='y', scilimits=(0,0), useMathText=False)
24 plt.show()

```

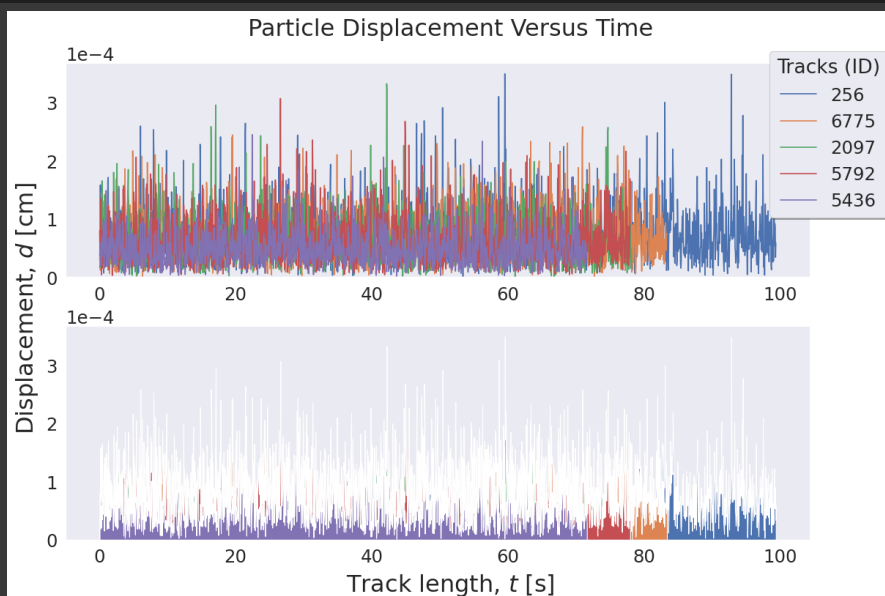

```
1 # units are wrong; need to be corrected in ImageJ before analysis.
2 # Good guide from JMU: https://www.jmu.edu/microscopy/resources/basic-image-processing-imagej.pdf
3 # Also read best practices for data analysis and presentation!
4
5 tracks_units
6 #file tracks.sort values(by=['TRACK DURATION'], ascending=False)
```


[illegible]

```
1 # diffusion coefficient
2
3 # changing data units
4 index = file_tracks.index.tolist()
5 df = pd.DataFrame(file_tracks.TRACK_ID)
6 df['TRACK_DURATION'] = file_tracks.TRACK_DURATION    #*0.05 # has been calibrated # 50 ms/frame = 0.05 s/frame (20.0 frames/s)
7 df['TRACK_DISPLACEMENT'] = file_tracks.TRACK_DISPLACEMENT*6.25    # 512px / 81.92e-6 m = 6.25 px/μm; 0.16 μm/px (?)
8 df['TRACK_DISPLACEMENT'] = df['TRACK_DISPLACEMENT']*1e-4    # 1 μm = 1e-4 cm
9
10 # r-squared and D
11 df['r2'] = (df.TRACK_DISPLACEMENT ** 2)
12 df['D'] = ( (df.TRACK_DISPLACEMENT ** 2) / ( 4 * df.TRACK_DURATION) )
13
14 # interpreting results: https://www.comsol.com/multiphysics/diffusion-coefficient
15 # In an aqueous (water) solution, typical diffusion coefficients are in the range of 1e-10 to 1e-9 m2/s
16 df
```

	TRACK_ID	TRACK_DURATION	TRACK_DISPLACEMENT	r2	D
	3	0.0	1.30	0.000018	3.206829e-10
	4	1.0	1.05	0.000057	3.285161e-09
	5	2.0	0.85	0.000061	3.697803e-09
	6	3.0	0.40	0.000127	1.612989e-08
	7	4.0	1.65	0.000090	8.069883e-09
...
	39540	39537.0	0.05	0.000032	1.046220e-09
	39541	39538.0	0.05	0.000080	6.335039e-09
	39542	39539.0	0.05	0.000057	3.298006e-09
	39543	39540.0	0.05	0.000328	1.073826e-07
	39544	39541.0	0.05	0.000067	4.435224e-09

39542 rows × 5 columns

```
1 # histogram using Seaborn + matplotlib
2
3 plot = sns.displot(data=df, x="D", kind="hist", kde=False, bins = 75, aspect = 1.5, legend=True)
4 plot.figure.subplots_adjust(top=0.9);
5 plt.xlim(-0.1e-6, None)
6 plot.figure.suptitle("Histogram of Diffusion Constants");
7 plot.set(xlabel=r'Diffusion Coefficient, $D$ $\left[ \mathrm{cm}^2/\mathrm{s} \right]$', ylabel='Absolute Frequency', xlim=(None, 2.5e-6));
```

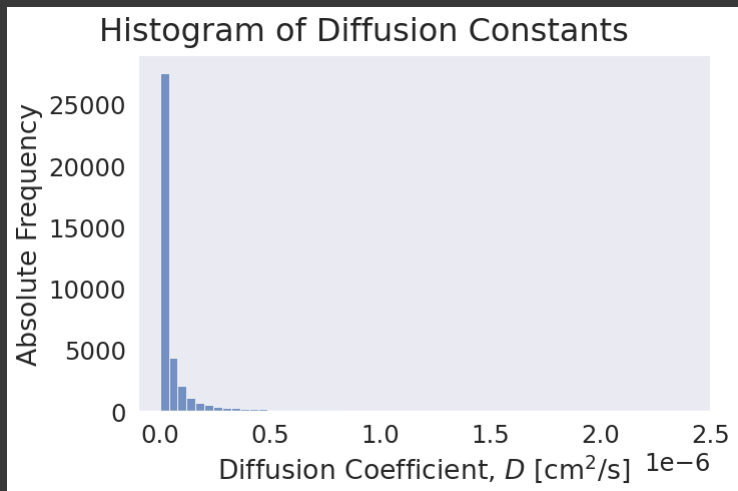

```
1 # histogram + kernel density estimate (KDE) plot
2
3 plot = sns.displot(data=df, x="D", kind="hist", kde=True, bins = 75, aspect = 1.5)
4 plot.figure.subplots_adjust(top=0.9);
5 plot.figure.suptitle("Histogram of Diffusion Constants With KDE Plot");
6 plt.xlim(-0.1e-6, None)
7 plot.set(xlabel=r'Diffusion Coefficient, $D$ $\left[ \mathrm{cm}^2/\mathrm{s} \right]$', ylabel='Absolute Frequency', xlim=(None, 2.5e-6));
8 plt.legend(labels=["KDE", "Histogram"]); # kernel density estimate (KDE) plot
```

```
1 # adding label to the df
2 df['Data_Series'] = 'pH 9 #04'
3
4 D_statistics = df.groupby(['Data_Series'])['D'].describe() #[['mean', 'std']]
5 D_statistics
```

	count	mean	std	min	25%	50%	75%	max
Data_Series								
pH 9 #04	39542.0	9.050378e-08	2.477131e-07	0.0	2.620137e-09	1.341095e-08	5.610246e-08	0.000003

▸ Michaelis-Menten Kinetics

pH 9.0

[] | 16 cells hidden

▸ Michaelis-Menten Kinetics

pH 7.8

[] | 5 cells hidden

▸ Code Snippets From Workshop

[] | 6 cells hidden

▸ Out of Scope

[] | 21 cells hidden